

Diagnostic Sodium Citrate Solution

診断用クエン酸ナトリウム液

Diagnostic Sodium Citrate Solution contains not less than 3.3 w/v% and not more than 4.3 w/v% of sodium citrate ($C_6H_5Na_3O_7 \cdot 2H_2O$: 294.10). The requirements as described for aqueous injections under Injections are applicable.

Method of preparation

Sodium Citrate	38 g
Water for Injection	a sufficient quantity
To make 1000 mL	

Prepare as directed under Injections, with the above ingredients.

No preservative may be added.

Description Diagnostic Sodium Citrate Solution is a clear, colorless liquid.

Identification Diagnostic Sodium Citrate Solution responds to the Qualitative Tests for sodium salt and citrate.

pH 7.0 – 8.5


Assay Pipet 5 mL of Diagnostic Sodium Citrate Solution, evaporate on a water bath to dryness, dry the residue at 180°C for 2 hours, and dissolve in 30 mL of acetic acid (100) by warming. Cool, and titrate with 0.1 mol/L perchloric acid VS (indicator: 3 drops of crystal violet TS). Perform a blank determination, and make any necessary correction.

$$\begin{aligned} \text{Each mL of 0.1 mol/L perchloric acid VS} \\ = 9.803 \text{ mg of } C_6H_5Na_3O_7 \cdot 2H_2O \end{aligned}$$

Containers and storage Containers—Hermetic containers.

Sodium Cromoglicate

クロモグリク酸ナトリウム

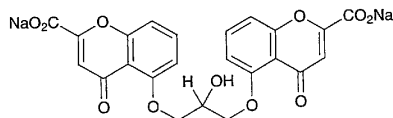

$C_{23}H_{14}Na_2O_{11}$: 512.33

Disodium 5,5'-(2-hydroxytrimethylenedioxy)bis(4-oxo-4H-1-benzopyran-2-carboxylate) [15826-37-6]

Sodium Cromoglicate contains not less than 98.0% of $C_{23}H_{14}Na_2O_{11}$, calculated on the dried basis.

Description Sodium Cromoglicate occurs as a white, crystalline powder. It is odorless and tasteless at first, and later develops a slightly bitter taste.

It is freely soluble in water, sparingly soluble in propylene glycol, very slightly soluble in ethanol (95), and practically insoluble in 2-propanol and in diethyl ether.

It is hygroscopic.

It gradually acquires a yellow color by light.

Identification (1) Dissolve 0.1 g of Sodium Cromoglicate in 2 mL of water, add 2 mL of sodium hydroxide TS, and boil for 1 minute: a yellow color is produced. After cooling, add 0.5 mL of concentrated diazobenzene sulfonic acid TS: a dark red color is produced.

(2) Determine the absorption spectrum of a solution of Sodium Cromoglicate in phosphate buffer solution, pH 7.4 (1 in 100,000) as directed under the Ultraviolet-visible Spectrophotometry, and compare the spectrum with the Reference Spectrum: both spectra exhibit similar intensities of absorption at the same wavelengths.

(3) Sodium Cromoglicate responds to the Qualitative Tests for sodium salt.

Purity (1) Clarity and color of solution—Dissolve 0.50 g of Sodium Cromoglicate in 10 mL of water: the solution is clear and colorless to pale yellow.

(2) Acid or alkali—Dissolve 2.0 g of Sodium Cromoglicate in 40 mL of freshly boiled and cooled water, add 6 drops of bromothymol blue TS, and use this solution as the sample solution. To 20 mL of the sample solution add 0.25 mL of 0.1 mol/L sodium hydroxide VS: a blue color is produced. To another 20 mL of the sample solution add 0.25 mL of 0.1 mol/L hydrochloric acid VS: a yellow color is produced.

(3) Heavy metals—Proceed with 1.0 g of Sodium Cromoglicate according to Method 2, and perform the test. Prepare the control solution with 2.0 mL of Standard Lead Solution (not more than 20 ppm).

(4) Oxalate—Dissolve 0.25 g of Sodium Cromoglicate in water to make exactly 50 mL, and use this solution as the sample solution. Separately, dissolve 0.049 g of oxalic acid dihydrate, exactly weighed, in water to make exactly 100 mL. Pipet 5 mL of this solution, add water to make exactly 100 mL, and use this solution as the standard solution. Pipet 20 mL each of the sample solution and the standard solution, add exactly 5 mL of iron salicylate TS to each solution, and add water to make 50 mL. Determine the absorbances of these solutions as directed under the Ultraviolet-visible Spectrophotometry using water as the blank: the absorbance of the sample solution at 480 nm is not smaller than that of the standard solution.

(5) Related substances—Dissolve 0.20 g of Sodium Cromoglicate in 10 mL of water, and use this solution as the sample solution. Pipet 1 mL of the sample solution, add water to make exactly 10 mL, pipet 1 mL of this solution, add water to make exactly 20 mL, and use this solution as the standard solution. Perform the test with these solutions as directed under the Thin-layer Chromatography. Spot 10 μ L each of the sample solution and the standard solution on a plate of silica gel with fluorescent indicator for thin-layer chromatography. Develop the plate with a mixture of methanol, chloroform and acetic acid (100) (9:9:2) to a distance of about 10 cm, and air-dry the plate. Examine under ultraviolet light (main wavelength: 254 nm): spots other than the principal spot from the sample solution is not more intense than the spot from the standard solution.

Loss on drying Not more than 10.0% (1 g, in vacuum, 105°C, 4 hours).

Assay Weigh accurately about 0.18 g of Sodium Cromoglicate, and dissolve in a mixture of 25 mL of propylene glycol

and 5 mL of 2-propanol by warming. After cooling, add 30 mL of 1,4-dioxane, and titrate with 0.1 mol/L perchloric acid-dioxane VS (potentiometric titration). Perform a blank determination, and make any necessary correction.

Each mL of 0.1 mol/L perchloric acid-dioxane VS
= 25.617 mg of C₂₃H₁₄Na₂O₁₁

Containers and storage Containers—Tight containers.
Storage—Light-resistant.

Sodium Iodide

ヨウ化ナトリウム

NaI: 149.89

Sodium Iodide, when dried, contains not less than 99.0% of NaI.

Description Sodium Iodide occurs as colorless crystals or a white, crystalline powder. It is odorless.

Sodium Iodide is very soluble in water, and freely soluble in glycerin and in ethanol (95).

It deliquesces in moist air.

Identification A solution of Sodium Iodide (1 in 20) responds to the Qualitative Tests for sodium salt and for iodide.

Purity (1) Clarity and color of solution—Dissolve 1.0 g of Sodium Iodide in 2 mL of water: the solution is clear and colorless.

(2) Alkali—Dissolve 1.0 g of Sodium Iodide in 10 mL of freshly boiled and cooled water, and add 1.0 mL of 0.005 mol/L sulfuric acid VS and 1 drop of phenolphthalein TS: no color is produced.

(3) Chloride, bromide and thiosulfate—Dissolve 0.20 g of Sodium Iodide in 5 mL of ammonia TS, add 15.0 mL of 0.1 mol/L silver nitrate VS, shake for a few minutes, and filter. To 10 mL of the filtrate add 15 mL of dilute nitric acid: no brown color appears. The solution has no more turbidity than the following control solution.

Control solution: To 0.30 mL of 0.01 mol/L hydrochloric acid VS add 2.5 mL of ammonia TS, 7.5 mL of 0.1 mol/L silver nitrate VS and 15 mL of dilute nitric acid.

(4) Nitrate, nitrite and ammonium—Place 1.0 g of Sodium Iodide in a 40-mL test tube, and add 5 mL of water, 5 mL of sodium hydroxide TS and 0.2 g of aluminum wire. Insert a pledget of absorbent cotton in the mouth of the test tube, and place a piece of moistened red litmus paper on the cotton. Heat the test tube on a water bath for 15 minutes: the evolved gas does not turn moistened red litmus paper to blue.

(5) Cyanide—Dissolve 0.5 g of Sodium Iodide in 10 mL of water. To 5 mL of this solution add 1 drop of iron (II) sulfate TS and 2 mL of sodium hydroxide TS, warm, and add 4 mL of hydrochloric acid: no green color develops.

(6) Iodate—Dissolve 0.5 g of Sodium Iodide in 10 mL of freshly boiled and cooled water, and add 2 drops of dilute sulfuric acid and 1 drop of starch TS: no blue color develops immediately.

(7) Heavy metals—Proceed with 2.0 g of Sodium Iodide

according to Method 1, and perform the test. Prepare the control solution with 2.0 mL of Standard Lead Solution (not more than 10 ppm).

(8) Barium—Dissolve 0.5 g of Sodium Iodide in 10 mL of water, add 1 mL of dilute sulfuric acid, and allow to stand for 5 minutes: no turbidity is produced.

(9) Potassium—Dissolve 1.0 g of Sodium Iodide in water, and add water to make 100 mL. To 4.0 mL of this solution add 1.0 mL of dilute acetic acid, shake, add 5.0 mL of a solution of sodium tetraphenylboron (1 in 30), immediately shake, and allow to stand for 10 minutes: the solution has no more turbidity than the following control solution.

Control solution: Dissolve 9.5 mg of potassium chloride in water, and add water to make 1000 mL. To 4.0 mL of this solution add 1.0 mL of dilute acetic acid, shake, and then proceed as directed above.

(10) Arsenic—Prepare the test solution with 0.40 g of Sodium Iodide according to Method 1, and perform the test using Apparatus B (not more than 5 ppm).

Loss on drying Not more than 5.0% (2 g, 120°C, 2 hours).

Assay Weigh accurately about 0.4 g of Sodium Iodide, previously dried, in an iodine flask, dissolve in 10 mL of water, add 35 mL of hydrochloric acid and 5 mL of chloroform, and titrate with 0.05 mol/L potassium iodate VS while shaking vigorously until the red-purple color of the chloroform layer disappears. The end point is attained when the red-purple color does not reappear in the chloroform layer within 5 minutes after the layer has been decolorized.

Each mL of 0.05 mol/L potassium iodate VS
= 14.989 mg of NaI

Containers and storage Containers—Tight containers.
Storage—Light-resistant.

Sodium Iodide (¹²³I) Capsules

ヨウ化ナトリウム (¹²³I) カプセル

Sodium Iodide (¹²³I) Capsules contain iodine-123 in the form of sodium iodide.

It conforms to the requirements of Sodium Iodide (¹²³I) Capsules in the Minimum Requirements for Radiopharmaceuticals.

Sodium Iodide (¹³¹I) Capsules

ヨウ化ナトリウム (¹³¹I) カプセル

Sodium Iodide (¹³¹I) Capsules are prepared by dispensing Sodium Iodide (¹³¹I) Solution into capsules and drying it.

Sodium Iodide (¹³¹I) Capsules conform to the requirements of Sodium Iodide (¹³¹I) Capsules in the Minimum Requirements for Radiopharmaceuticals.